BLACKROLL®

INFINITY SPARE BATTERIES 2.0

Instruction Manual

EN | DE | FR | NL | IT | PL | CN

Hi!

Congratulations on your new BLACKROLL® INFINITY SPARE BATTERIES 2.0.





Training and recovery are the key to your well-being and performance. With the INFINITY SPARE BATTERIES 2.0 you can use your COMPRESSION BOOTS even longer. In addition to the increased battery capacity, they have a function for continuous operation. Before using the product, please read the instructions for use to ensure correct use of the product,

Safety Instructions

- + Persons with reduced physical, sensory or mental capabilities and any person under the age of 18 must be supervised during use.
- + Only use the device with the charger supplied.
- + Keep the device away from heat sources, children and pets.
- + Protect the device from moisture, dust, aggressive liquids and vapors.
- + Should moisture nevertheless penetrate the device, switch it off immediately and, if necessary, disconnect the mains plug from the socket.
- + Do not remove any screws and do not attempt to dismantle the appliance.
- + Only pass the appliance on to third parties with the operating instructions.
- + Any warranty claim is void if the appliance is misused.

Instructions for use

- A** USB-C connection
- B** Battery indicator
- C** Output end of the cable
- D** Input end of the cable
- E** Power supply unit
- F** Charge status

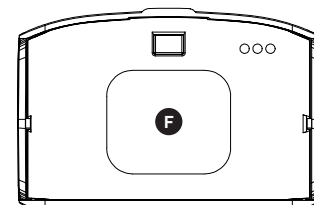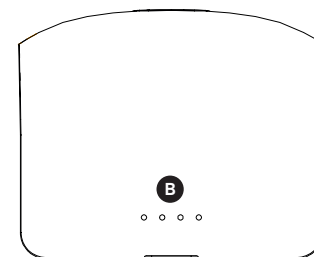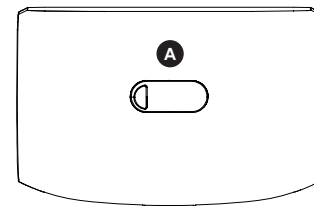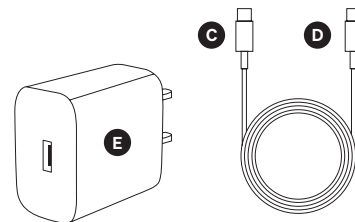Before first use

Charge the batteries fully before first use. Charging time approx. 3 hours. The battery is fully charged when all 4 LEDs on the battery status indicator (B) are permanently lit.

To do this, connect the output end of the cable (C) to the USB-C port (A) and the input end of the cable (D) to the supplied power supply unit (E).

Your batteries are empty as soon as the battery status indicator (B) flashes during operation. If you want to continue using your COMPRESSION BOOTS even though the batteries are empty, connect the supplied charging cable with mains adapter as described above. As soon as the status indicator (B) lights up continuously, the batteries are being charged. You can charge the batteries at any time during operation.

With fully charged batteries, the operating time is around 3.5 hours, depending on the mode and intensity of use.

If you want to check the battery status in between, you can do this by holding down the F button (charge status). Each lit LED corresponds to 25% battery charge. If the battery level is depleted to less than 10%, the first indicator light flashes and the battery must be recharged.

Cleaning

Do not clean the battery while it is charging. Do not use water to clean the device, only a damp cloth. Do not use corrosive cleaning agents or solvents. Textile parts are not suitable for cleaning in the washing machine.

Technical specifications

Battery:	Rechargeable and removable 2000mAh lithium-ion battery
Work in watt hours:	14.8
Battery capacity in mAh:	2000
Maximum charging current in A:	2
Maximum charging voltage in V:	5
Can the battery be replaced:	Yes
Type of charging connection:	USB-C
Nominal voltage:	DC7.4V
Is a charging cable supplied?	Yes
Power adapter included?	Yes, use a 5V 4A power adapter (Note: This specific adapter must be used)
Protection class:	IP-20
Battery life:	Up to 3.5 h, depending on mode and usage intensity
Battery charging time:	3h

Disposal

The unit contains electronic components and a battery. It must not be disposed of in normal household waste. All electrical and electronic equipment must be disposed of at government designated sites.

Warranty

The seller grants the buyer, if he is a consumer, the statutory warranty in accordance with the legal requirements for a period of 2 years from the purchase of the item. The Seller's warranty obligations shall initially be limited to subsequent performance, whereby the Seller shall be free to choose whether to deliver a defect-free item or to repair the defective item. If subsequent performance is unsuccessful, the buyer shall be entitled to the further statutory claims. The seller rejects warranty claims for damage caused by accident, obviously incorrect use, misuse or improper transport; which occurred as a result of repairs carried out by unauthorised personnel; which occurred as a result of use and handling not in accordance with the operating instructions; which occurred as a result of improper handling of the purchased item. Furthermore, claims for damages by the purchaser are excluded. Excluded from this are claims for damages by the buyer arising from injury to life, limb or health as well as liability for other damages based on an intentional or grossly negligent breach of duty by the seller, his legal representatives or vicarious agents. Also unaffected is the liability for breaches of duty, the fulfilment of which makes the proper execution of the contract possible in the first place and on the observance of which the buyer may regularly rely. In the event of a slightly negligent breach of such contractual obligations, the Seller shall only be liable for the foreseeable damage typical for the contract, unless it is a matter of claims for damages by the Customer arising from injury to life, limb or health.

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BLACKROLL AG
WSVK Oederan GmbH